



Helicobacter pylori quick test

Effortless biopsy testing



Ordering details

REF	Product
602015	Helicobacter pylori quick test, 50 tests
602017	BIOHIT Helicobacter pylori Control, 1 bottle



Contact our team for more information.



Biohit Oyj
Laippatie 1, 00880 Helsinki, Finland
Tel. +358 9 773 861
info@biohit.fi

www.biohithealthcare.com

490062en-4.0 ©BIOHIT OYJ 09/2026



BIOHIT
Innovating for Health

Easy testing from a biopsy specimen

The *H. pylori* quick test is a one-step test method to detect *H. pylori* infection from a biopsy sample during gastroscopy. The *H. pylori* quick test can be used to diagnose *H. pylori* infection or to determine the success of eradication therapy.

- Testing and reporting during gastroscopy
- One-step test procedure
- Positive results in 2 min (neg. 30 min)
- Sensitivity 100%, Specificity 93%, Accuracy 96% (ref. to histology)¹
- Clear indicator color – easy visual interpretation


Using the *H. pylori* quick test is an easy one-step procedure. The biopsy specimen is immersed into the gel medium, and if *H. pylori* urease is present in the specimen, a red color develops in the gel. Positive results are ready in only 2 minutes, and the final confirmation of a negative result can be made at 30 minutes. Interpretation of the indicator color is simple and does not require any specialist training.

Product specifications

Sample type	Gastric biopsy specimen
Test principle	<i>H. pylori</i> urease activity
Performance	Sensitivity 100%, Specificity 93% (ref. to histology)
Storage	2 to 8°C
Test results	2-30 minutes
Quality control	BIOHIT <i>Helicobacter pylori</i> Control available
Kit components	Test plate with gel medium, Instructions for Use


One-step test procedure

1. Place the biopsy specimen into the gel and incubate for 2 minutes.

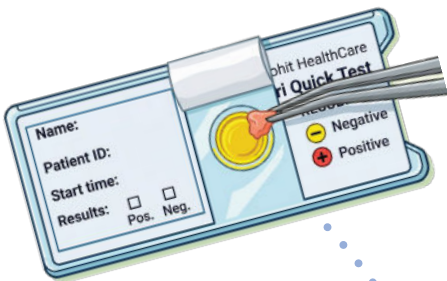


2. Read the test results:
2 minutes for positive results

RESULTS IN

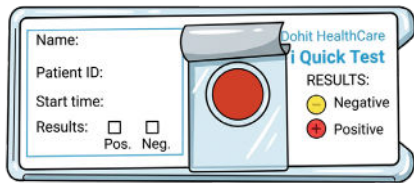

30 minutes for negative results

RESULTS IN

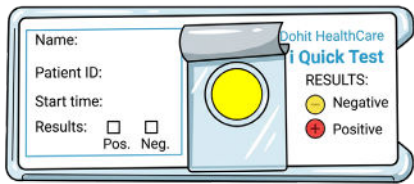

Reference

1. Vauhkonen M, et al. 2008. A new simplified urea test for *H. pylori* urease activity for gastric biopsy. *Helicobacter*. 13 (5): 473.