

BIOHIT TOTAL 25OH VITAMIN D

ELISA kit



Ordering details:

Cat. No.	Item
602310	BIOHIT Total 25OH Vitamin D ELISA, 96 tests

Biohit Oyj

Laippatie 1
00880 Helsinki
Finland
Tel. +358 9 773 861
Fax +358 9 773 86 200
info@biohit.fi

Biohit HealthCare Srl

Via Boncompagni, 3
20139 Milano
Italy
Tel +39 02 38238113
Fax +39 02 38236521
info.italy@biohit.fi

Biohit HealthCare Ltd.

Pioneer House, North Rd
Ellesmere Port, CH65 1AD,
United Kingdom
Tel. +44 151 550 4 550
Fax. +44 151 550 4 551
info@biohithealthcare.co.uk

www.biohithealthcare.com

BIOHIT HealthCare
Innovating for Health

BIOHIT HealthCare
Innovating for Health



Product specifications

- Detects both 25OH vitamin D2 and D3 for a clinically meaningful assessment of vitamin D status
- Calibrated against the ID-LC/MS-MS Reference Measurement Procedure
- User-friendly and fully automatable assay protocol



In the USA and Japan for research use only. Not available in Iran and Spain.

490097EN 10/2015

BIOHIT TOTAL 25OH VITAMIN D

ELISA kit

FOR A CONCLUSIVE DETERMINATION OF
VITAMIN D STATUS

BIOHIT TOTAL 25OH VITAMIN D

ELISA kit

Multifunctional vitamin D

Vitamin D has multiple roles in the human body. In addition to its well established role in the regulation of calcium absorption and promoting bone growth, it is recognized for other health benefits including reducing risk of diseases such as type 1 diabetes and common cancers.

Clinical background

Two major forms of vitamin D are important in the human body, vitamin D2 (ergocalciferol) and vitamin D3 (cholecalciferol). Vitamin D3 is the form synthesized in the skin in response to sunlight UVB exposure yet can also be obtained from animal based foods. Vitamin D2, on the other hand, is the form synthesized by plants and is obtained from plant derived foodstuff. Both are metabolized in the liver to their respective 25OH vitamin D3 and D2 forms which are in turn converted to their active form in the kidneys.

The best indicator of vitamin D status is the serum concentration of 25OH vitamin D. For a correct diagnosis of vitamin D deficiency, insufficiency or intoxication, the assay must recognize both the D2 and D3 forms. As both forms have clinical relevance, the total concentration provides the information that can be acted upon.

BIOHIT Total 25OH D ELISA kit


The BIOHIT Total 25OH vitamin D ELISA kit is a quantitative immunoenzymatic assay. While detecting both 25OH vitamin D2 and D3, the kit provides clinically relevant information on the vitamin D status. Reliability of the results is ensured by validation against the ID-LC/MS-MS Reference Measurement Procedure (Ghent method) as approved by the Vitamin D Standardization Program (VDSP) with $R > 0.97$.



The assay protocol


The assay protocol includes just six steps making it easy to perform and reducing the hands-on time. Featuring pretreatment technology performed in the assay microtiter plate well, any open ELISA platform can be easily programmed to run the assay from the start.

1. Add 50 μ l of sample, calibrators and controls and 150 μ l of incubation buffer to each well

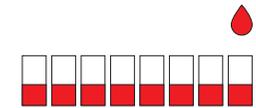


Incubate 120 min at RT with shaking
→ Wash plate

2. Add 200 μ l conjugate

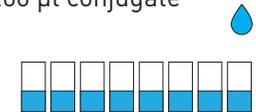


Incubate 30 min at RT with shaking
→ Wash plate

3. Add 100 μ l TMB

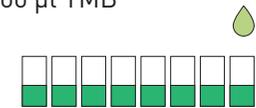

4. Incubate 15 min at RT with shaking

5. Add 100 μ l STOP

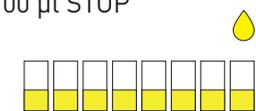

6. Read Absorbance at 450 nm