



#### Product specifications:

**Sample:** gastric biopsy specimen

**Test principle:** *H. pylori* urease activity

#### Kit components:

- Test tube/plate
- Reagent
- Instructions for Use

**Storage:** 18 to 25°C

**Quality control:** BIOHIT *Helicobacter pylori*  
Control available



Not available in the United States. In Japan for research use only.  
Check availability for other markets.

#### References

1. Vaira D, Vakil N, Gatta L, Ricci C, Perna F, Saracino I, Fiorini G, Holton J. Accuracy of a new ultrafast rapid urease test to diagnose *Helicobacter pylori* infection in 1000 consecutive dyspeptic patients. *Aliment Pharmacol Ther* 31 (2010), 331-338.

#### Ordering details:

##### BIOHIT *Helicobacter pylori* UFT300 with plate

| REF | Product |
|-----|---------|
|-----|---------|

|           |                                                   |
|-----------|---------------------------------------------------|
| 602005PLA | BIOHIT <i>Helicobacter pylori</i> UFT300, 5 tests |
|-----------|---------------------------------------------------|

|           |                                                    |
|-----------|----------------------------------------------------|
| 602019PLA | BIOHIT <i>Helicobacter pylori</i> UFT300, 50 tests |
|-----------|----------------------------------------------------|

##### BIOHIT *Helicobacter pylori* UFT300 with tube

| REF | Product |
|-----|---------|
|-----|---------|

|        |                                                    |
|--------|----------------------------------------------------|
| 602019 | BIOHIT <i>Helicobacter pylori</i> UFT300, 50 tests |
|--------|----------------------------------------------------|

|        |                                                     |
|--------|-----------------------------------------------------|
| 602021 | BIOHIT <i>Helicobacter pylori</i> UFT300, 100 tests |
|--------|-----------------------------------------------------|

| REF | Product |
|-----|---------|
|-----|---------|

|        |                                           |
|--------|-------------------------------------------|
| 602017 | BIOHIT <i>Helicobacter pylori</i> Control |
|--------|-------------------------------------------|

Contact our team for more information.



Biohit Oyj  
Laippatie 1, 00880 Helsinki, Finland  
Tel. +358 9 773 861  
info@biohit.fi

[www.biohithealthcare.com](http://www.biohithealthcare.com)

## BIOHIT *Helicobacter pylori* UFT300

Ultra-fast and effortless biopsy testing




490093en-05.2 ©BIOHIT OYJ 08/2025


**BIOHIT**  
Innovating for Health

## Simple test procedure with plate


1. Place 1 biopsy into the well of the plate.

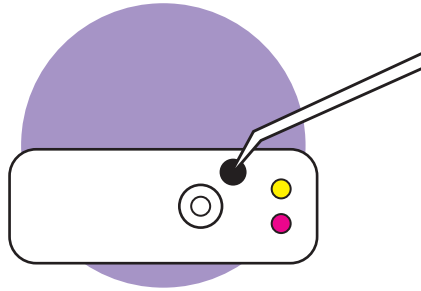

2. Add one drop of the reagent solution. Cover the plate with label and shake for 5 seconds.

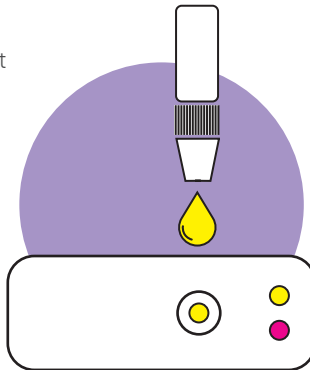

3. Read after 5 minutes.

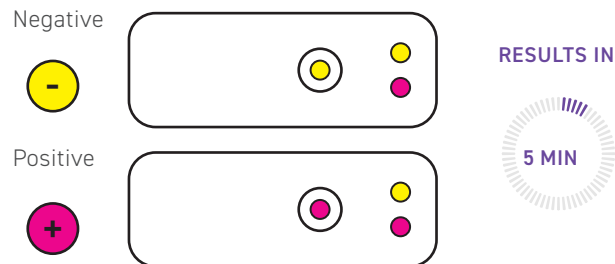

## Ultra-fast *H. pylori* detection from a biopsy

The BIOHIT *H. pylori* UFT300 is a true quick test for the detection of *H. pylori* from a biopsy. The biopsy taken during gastroscopy can be tested immediately to diagnose *H. pylori* infection or to check the success of eradication therapy. The result is ready in just 5 minutes enabling diagnosis and reporting in endoscopy. The UFT300 has excellent sensitivity and specificity which makes it a highly reliable tool for diagnostics.


- Ready-to-use test kit
- Results ready in 5 min
- Storage at room temperature
- Sensitivity 94,5%, Specificity 100%<sup>1</sup>
- Test and report immediately

## Biopsy testing could not be easier


The BIOHIT *H. pylori* UFT300 is easy and effortless to use. The biopsy specimen is placed into the test tube/well of the plate and mixed with the test reagent. Color change to magenta indicates the presence of *H. pylori* in the specimen. The procedure is safe and hygienic for the user, and the interpretation of the result doesn't require specialist training.

## Simple test procedure with tube


1. Place 1 biopsy into the tube.

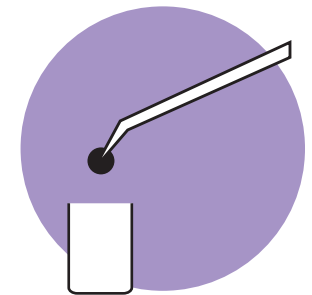

2. Add one drop of the reagent solution. Cap the tube and shake for 5 seconds.

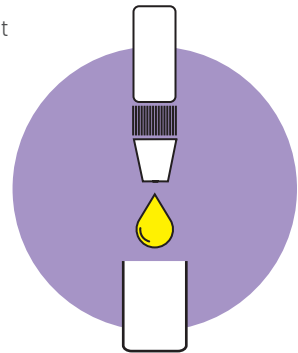

3. Read after 5 minutes.

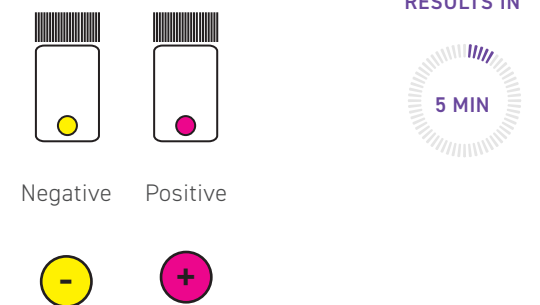

"The result is ready in just 5 minutes."