



BIOHIT Celiac quick test



RESULTS IN



Product specifications:

- Sample: fingerprick blood, plasma or serum
- Test principle: immunological detection of anti-tTG antibodies (IgA, IgG and IgM)
- Kit components: lancets, micropipettes, test cassettes, diluent buffer, Instructions for Use
- Storage: 2 to 8°C



Ordering details:

REF	Product
602070	BIOHIT Celiac quick test, 20 tests
602080	BIOHIT Celiac quick test, 50 tests

Contact our team for more information.



Biohit Oyj
Laippatie 1, 00880 Helsinki, Finland
Tel. +358 9 773 861
info@biohit.fi

www.biohithealthcare.com

BIOHIT Celiac quick test

Reliable first-line diagnostics



Not available in the United States. In Japan for research use only.
Check availability for other markets.

490068en-05.2 ©BIOHIT OYJ 08/2025

BIOHIT
Innovating for Health

Accurate results from fingerprick blood

BIOHIT Celiac quick test enables accurate and patient-friendly testing of celiac disease from only a drop of blood. The test is based on the detection of anti-tissue transglutaminase antibodies (tTG); (IgA, IgG, IgM), which delivers superior sensitivity, especially in patients with IgA deficiency.

- Results in 10 minutes
- Fingerprick sample collection
- Easy-to-use lateral flow quick test



Easy testing at the Point-of-Care

The blood sample for BIOHIT Celiac quick test is taken as a fingerprick sample. The sample is placed on the test cassette along with a buffer solution. The test result is clearly visible within 10 minutes.


“Fast results in any clinical setting.”

Easy sampling - fast results

1. Collect the fingerprick sample.

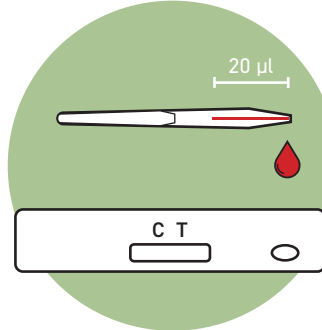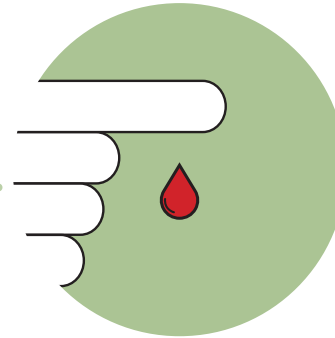

2. Dispense blood sample (20 µl) and two drops of buffer on the test cassette.

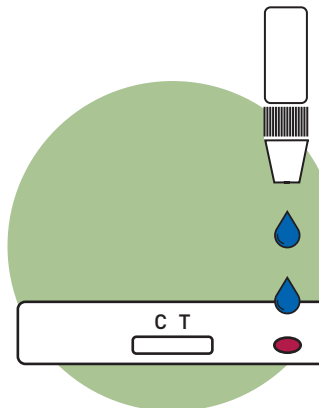

3. Read after 10 minutes

RESULTS IN



Convenient, even
for pediatric patients

NEGATIVE



POSITIVE

