



BIOHIT ColonView® quick test (FIT)

For early diagnosis of colorectal cancer



Product specification

Sample type	Stool
Test principle	Immunological detection of Hb and Hb/Hp complex
Performance	Adenoma and carcinoma sensitivity 97% and specificity 85%
Storage	2 to 30°C
Test results	15 minutes
Quality control	BIOHIT ColonView® QT Control available
Kit components	Test cassettes, sample collection package (sample collection tube, specimen collection paper, and sample collection instruction) and Instructions for Use

Reference

1. Vasilyev S, Smirnova E, Popov D et al. Anticancer Res 2015; 35:2873-2880.

Ordering details

REF	Product
602250.02	BIOHIT ColonView® quick test, 30 tests
602390	BIOHIT ColonView® QT Control



Contact our team for more information.



Biohit Oyj
Laippatie 1, 00880 Helsinki, Finland
Tel. +358 9 773 861
info@biohit.fi

www.biohithealthcare.com

490089en-4.0 ©BIOHIT OYJ 09/2026

BIOHIT
Innovating for Health

Detects bleeding derived from lower GI tract

Detection of occult blood in feces is essential in the early diagnosis of colorectal cancer. The BIOHIT ColonView® quick test is a qualitative *in vitro* immunochromatographic test aiding the diagnosis of lower gastrointestinal disorders, such as colorectal cancers and large adenomas that bleed, by detecting human hemoglobin (Hb) and hemoglobin/haptoglobin complex (Hb/Hp) in stool samples. With this combination of two biomarkers, the ColonView test can be used for detecting bleeding originating in the proximal and distal colon. This makes ColonView a superior tool to detect bleeding cancer tumors and their precursors.

- Detects bleeding from proximal and distal lesions
- Adenoma & carcinoma: sensitivity 97% and specificity 85%¹
- Hygienic sampling, transportation and analysis
- No dietary restrictions before sample collection
- Results in 15 min

One sample for dual testing

The innovative design of the BIOHIT FAEX® Sample system eliminates patient-errors and ensures that there will always be the correct sample amount. The sampling tube contains all the necessary reagents and no additional accessories are needed for the analysis. One sample is enough for the testing of both Hb and Hb/Hp complex, and the result is ready in 15 minutes.

Simple test procedure - fast results


1. Collect the sample.



RESULTS IN



3. Read the test results after 15 minutes.

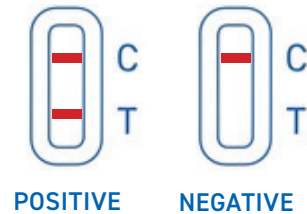

2. Add diluted sample to the test cassette.

